


# Ozonobaric P

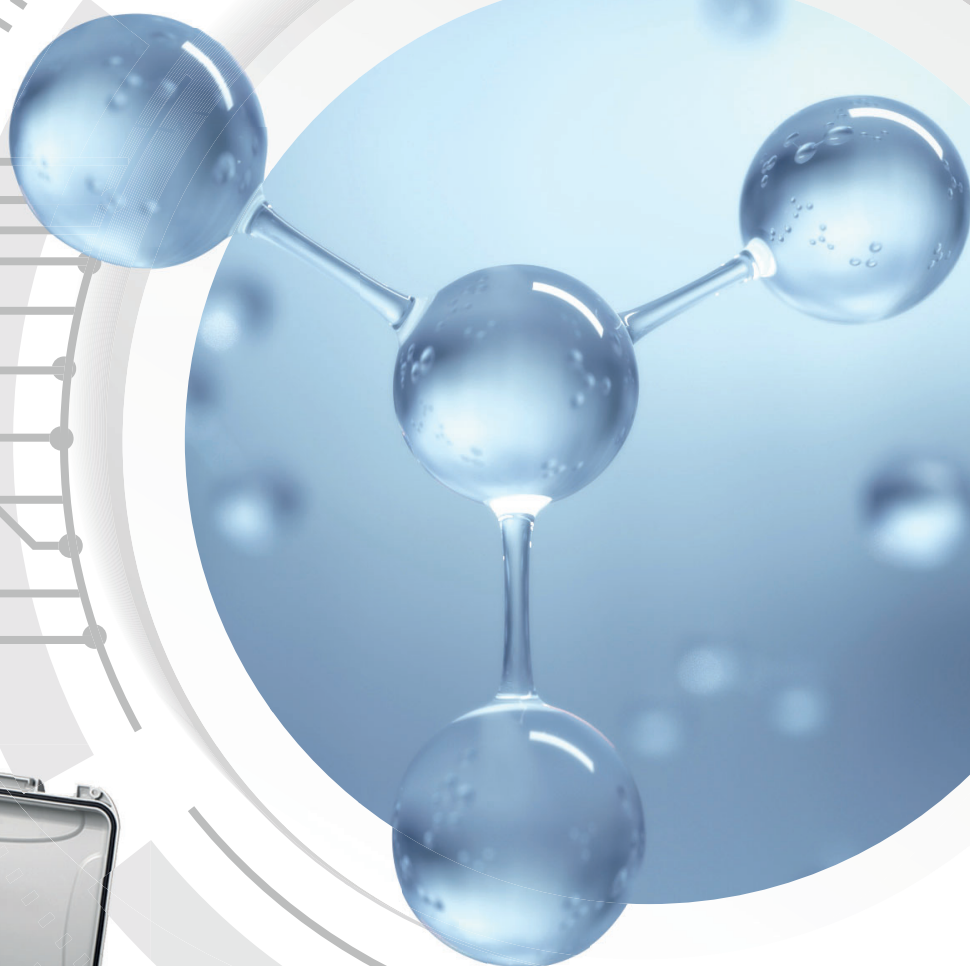

# Technical features

**Concentration range** 1 µg/Nml to 80 µg/Nml Standardized measurement units.

**Output** 3 g/h.

**Continuous concentration adjustment** in steps of 1 µg/Nml.

**Concentration measurements**, dual, photometric and indirect (ACM). Maintenance-free.

**Precision of measurements** accuracy in concentration.

**Adjustable flow rate** in manual and dynamic modes from 10 l/h to 50 l/h.

**Modes of operation** syringe, dose, bag, continuous, insufflation, vacuum.

**Integrated vacuum pump**, Vacuum up to -700 mbar (-10.15 PSI).

**Independent ozone and vacuum inlets**, standard Luer-lock system.

**Contactless capacitive keypad and LCD display.**

**Power supply** 100-230 VAC/ 50 - 60 Hz/ 145 VA.

**Oxygen inlet by quick tapping**. Inlet pressure from 1.5 bar to 4.5 bar.

**Weight** 9 Kg - 19.8 Lb (without O<sub>2</sub> bottle).

**Dimensions** 50.3 cm x 40.6 cm x 19.3 cm (19.8 x 16 x 7.6 Inch).

**Certifications**, CE, ( MD II B class ), CSA, RoHS

The **OZONOBARIC P** medical ozone generator is a portable system, integrated in a **IP67 WATERPROOF** case, **NK-7 IMPACT RESISTANT**, valve with automatic pressure adjustment and PowerClaw closure, ideal for Ozone applications in Human and Veterinary Medicine.

Suitable for any type of ozone therapy:

- Major Auto Hemotherapy (MAHT) Minor Auto Hemotherapy (mAHT) intraarticular and intramuscular infiltrations, new specific mode for vaginal and rectal insufflation, topical treatments in automatic (open and closed circuit) and liquids to be ozonized. Designed with the most rigorous standards of the electro-medical industry, incorporating advanced technology with a double measurement system: direct photometric system by LED (complies with RoHS regulations) to continuously measure the concentration produced and ACM indirect high precision system.

The double circuit of internal ozone catalysts for complete isolation between generator and patient, avoids possible cross-contamination.

The compact and high precision tube, inert to the ozone action, and high mechanical strength is manufactured with state-of-the-art materials and consist of a high voltage generator™ by dielectric discharge.

Our high voltage generator has been specifically designed, taking advantage of the latest technological innovations from I+D+I of our XR generators. It also has an integrated vacuum pump for independent use.

High efficiency with minimum O<sub>2</sub> consumption, free of annual low maintenance cost.



**Easy to carry**



**Safety and precision**

**SEDECAL S.A.**

C/ Pelaya, 9-13  
28110 Algete,  
Madrid (SPAIN)  
+34 916 280 544  
O3@sedecal.com  
www.o3.sedecal.com



CAT\_007\_OZO\_23\_EN\_OZONOBARIC P\_REV#1

**SEDECAL USA**

7555 Caldwell Avenue  
Niles, IL. 60714, USA  
+847.394.6960  
O3@sedecal.com  
www.o3.sedecal.com